

In step with Capintec

The Newsletter of Capintec, Inc.

September/October 2003

Stop By & Visit Us:
Booth # 336 at the ASNC 8th
Annual Meeting in Indianapolis

contents

Generation II Cardiac Vest®

page 2

Nuclear Cardiology Hot Lab Packages

page 3

CRC®-15 Enhanced Software

page 4

Visit Capintec at the ASNC Annual Meeting in Indianapolis

The American Society of Nuclear Cardiology will be holding its 8th Annual Symposium and Scientific Session (September 11-14, 2003) at the Indiana Convention Center. Capintec will be there to support the Society and to exhibit its full line of products designed for Nuclear Cardiology.

Attendees are welcome to take this opportunity to stop by, look and learn about our exciting products and see why Capintec is at the heart of Nuclear Cardiology.

Be sure to visit us in **Booth # 336** in order to see Capintec's products and pick up your free copy of our new Radiology and Nuclear Medicine catalogs

Capintec.com "Enhanced"

Click on www.capintec.com and visit our ENHANCED website and view our full line of quality Nuclear Medicine, PET/511, Radiation Therapy & Diagnostic Radiology products. If you're looking for more information or need a quote, try our convenient shopping "basket" feature. Also, find our meeting schedules, special offers & featured products, all with the click of your mouse.



Home

About Us

Capintec News

Products

Request Info

Support

Contacts

Capintec, Inc.



Meeting Schedule


Featured Product

Special Offers

Capintec, Inc.
6 Arrow Road, Ramsey, NJ 07446
Phone: (201) 825-9500 / (800) 631-3826
Fax: (201) 825-4829

Capintec... One Stop for your Nuclear Cardiology Hot Lab Solutions

If you are setting up a new department or office, remodeling an existing hot lab or seeking a product to meet a specific requirement, Capintec can help. We can customize products to your specifications and create the hot lab package that meets your needs.

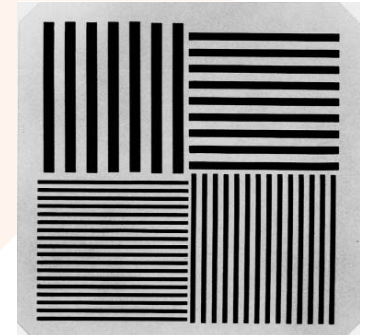

To request information and a quotation, call 1-800-631-3826 or email us at getinfo@capintec.com

CRC[®]-15 Enhanced Software

Add more functionality to your Capintec CRC[®]-15R Dose Calibrator with CRC[®]-15R ENHANCED SOFTWARE, all of the same great features of the CRC[®]-15R plus a few new ones.

The CRC[®]-15R Dose Calibrator has a long history of accuracy and reliability, so take that proven track record to the next level. The CRC[®]-15R ENHANCED SOFTWARE will still perform the standard Daily QC protocol, including zero adjust, background subtraction, system test, accuracy and constancy. But, it can now calculate dose tables, perform geometry tests (for between 2 and 10 measurements), do radiochromatography, calculate linearity (using standard, Lineator or Calicheck methods) and maintain a radiopharmaceutical inventory of kits or unit doses, with all acquired data being printable.

The Enhanced Software will solidify the reputation of the CRC[®]-15R as the workhorse of the Nuclear Medicine community, on the job 24/7/365.

The CRC[®]-15R Enhanced Software is an available option in all new CRC[®]-15R Dose Calibrators or it can be added to your existing CRC[®]-15R.

- Inventory control
- Dose table calculations
- Automated linearity protocol (standard, Lineator or Cali-check methods)
- Built in geometry testing (between 2 and 10 measurements)
- Radiochromatography QC package

STANDARD LINEARITY

#	MEASURED	%VAR
1	50.4mCi	2.10
2	36.5mCi	1.85
Any Key to Continue		

Standard Linearity
Test Results Screen

INVENTORY

1. ADD
2. KIT
3. WITHDRAW
4. PRINT
5. DELETE

Inventory Menu